


**dx** dental x the sterilization company

Committed to quality and innovation

STERILINE 1023

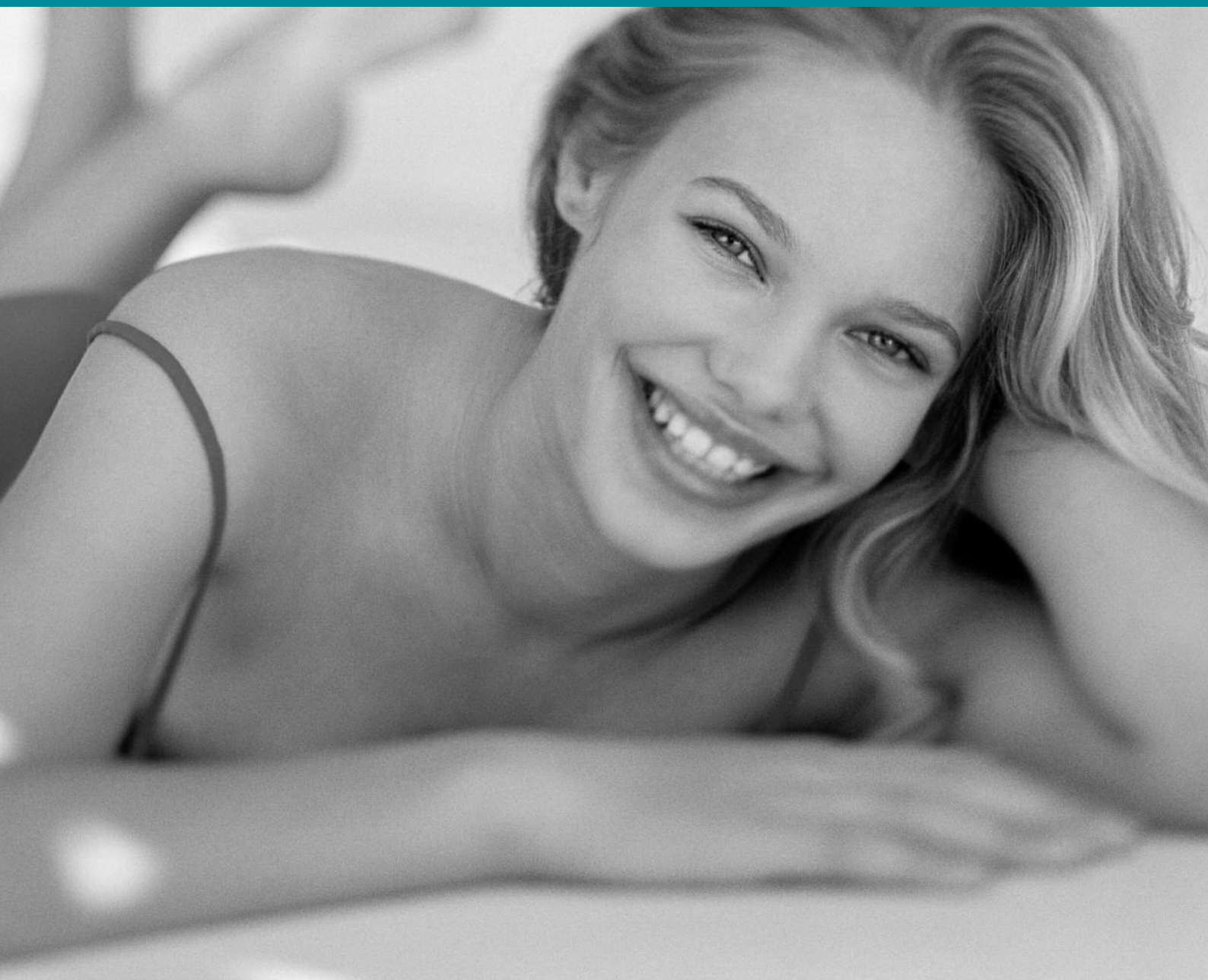

# STERILINE 1023

## Sealing system with adjustable temperature control

The Steriline 1023 system enables surgical instruments and handpieces to be wrapped before being put in the autoclave so as to maintain sterility after sterilization has been completed.

The 10-12 mm seal operated by the Steriline 1023 further accomplish the seal thickness requested by the norm EN ISO 868-5 and it grant the higher safety standards.

A potentiometer allows to adjust the sealing temperature in order to grant the perfect seal of every type of paper. The Steriline 1023 sealing technology grants the higher reliability and an easy use to the staff.

A visual indication and a sound signal will inform the user that sealing have been finished as requested by the norm EN ISO 11607-2. Such feature prevents from the element overheating.

The choice of advanced components and materials and the strict quality manufacturing and testing controls ensure the highest reliability and an easy operation to the clinical staff. The compact design with smooth and rounded surfaces make the unit easier to clean.


- Automatic 10-12 mm seal
- Digital temperature control
- Scorch protection
- Blade snapback mechanism after cutting
- Reel support
- Committed to quality and innovation

The instruments preparation has never been so easy




## the safe sealing in compliance with the norm

- Automatic start of sealing as soon as the lever is depressed, with visual display of operation.
- The temperature control permits to adjust the temperature in order to seal different type of paper.
- The snapback blade permits to size every package and reduce paper waste.

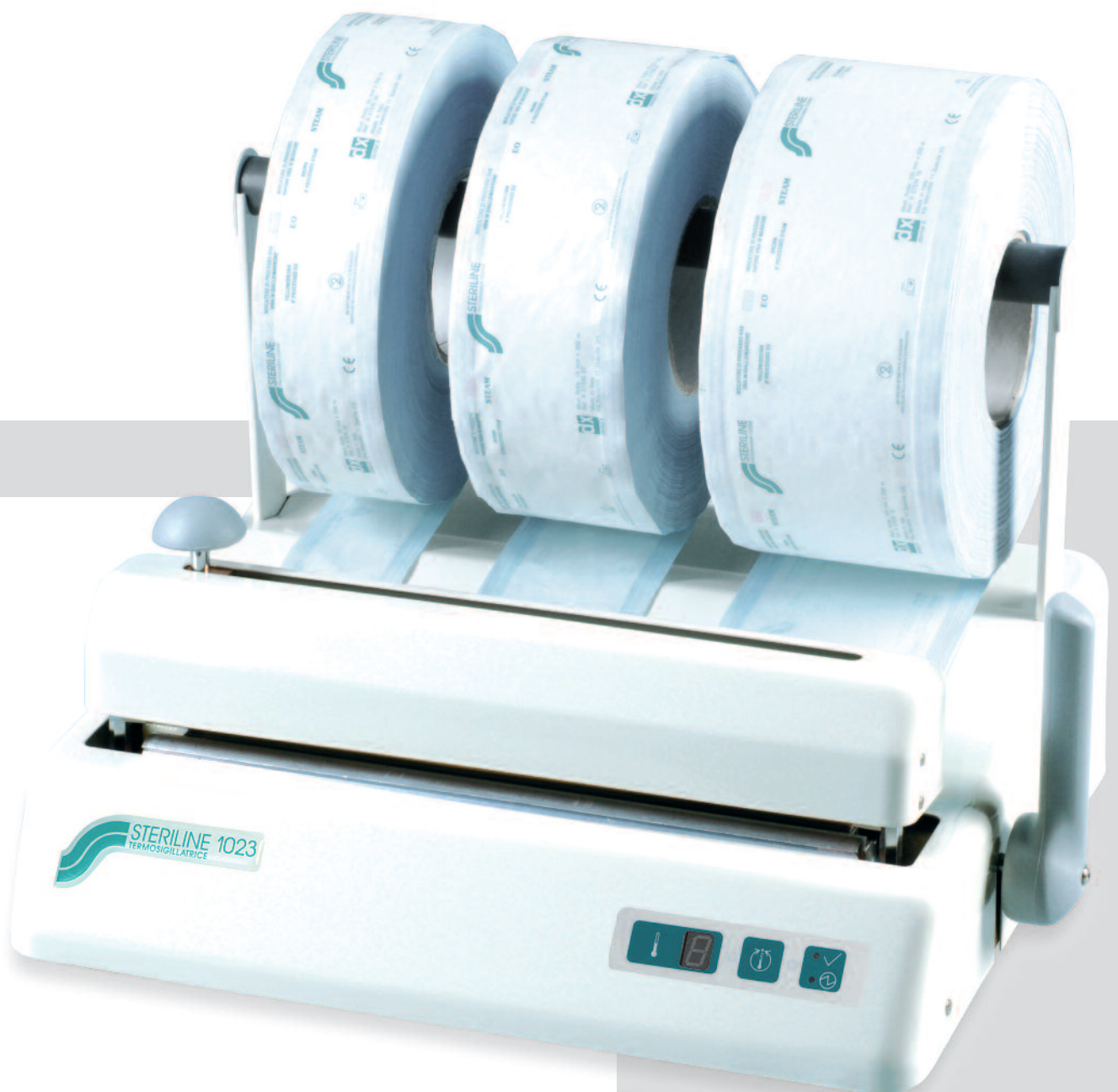

Technical features	Steriline 1023
Seal	10-12 mm
Blade snapback mechanism	■
Automatic start of sealing as soon as the lever is depressed	■
Visual and sound alarm at the end of the seal	■
Reel support	■
Lever scorch protection	■
Metallic bodywork	RAL 9003
Dimension (HxDxW)	335x323x421
Weight	8,5 kg
Electric Power	230V 50/60Hz
Consumption Max	600W
Wall support	■
Pre-heating time	30 sec
Heating system	Adjustable
Visual indication	digit
CE	■
ISO 11607-2	■
ISO 868-5	■

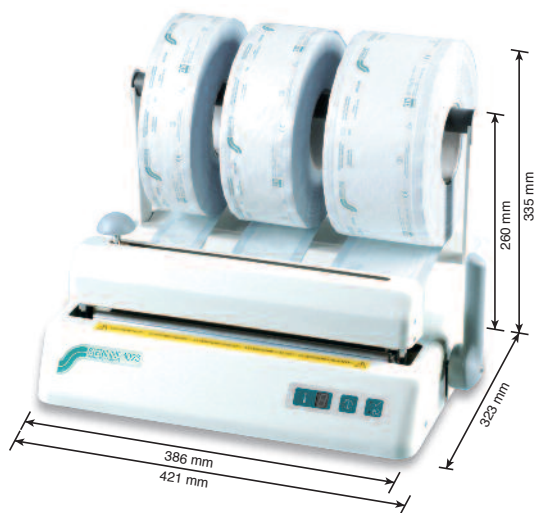

#### STERILINE MED-V

This new rotary sealer permits to operate the instrument packaging automatically and rapidly.

Steriline Med-V is supplied with the validation in compliance with the norm EN 11607-2.



Dental X reserves the right to make alterations without prior notice.

Tipo-Ito ABALTI G.F. - Creazzo (VI) • Code: GB\_STERILINE1023\_03-11

[www.dentalx.it](http://www.dentalx.it)



dental X spa  
 via marzotto 11  
 36031 dueville vicenza  
 tel. +39 0444 367400  
 fax +39 0444 367436  
 e mail [dentalx@dentalx.it](mailto:dentalx@dentalx.it)  
 internet [www.dentalx.it](http://www.dentalx.it)